

# XFamily<sup>2</sup>

**Generation<sup>2</sup> - Great functionality and look**



# You can trust on Breezy!

The new, improved Breezy XFamily<sup>2</sup> will give your business a boost.  
Efficiency which pays off.

Breezy is pure reliability.  
Robust materials, safe travel on any terrain, an excellent  
price-performance ratio.  
Improved functionality for an improved positioning of the  
user and more safety.  
The whole XFamily<sup>2</sup> is crash tested according ISO 7176-19.




# The Breezy XFamily<sup>2</sup>!

Optimise your fleet management.  
 Nine wheelchairs – one concept. In next to no time, you can swap parts and carry out individual adjustments.

## The advantages for you at a glance:

- Sophisticated specialist retail range – you will be able to offer the full spectrum of products in the field of re-use
- Common parts ensures excellent profitability
- Easy to re-use offers a high flexibility
- Easy re-configuration, adjustment and retro-fitting – it couldn't be any quicker or easier
- Low stock levels – minimal capital tied up
- Less handling effort – quicker results

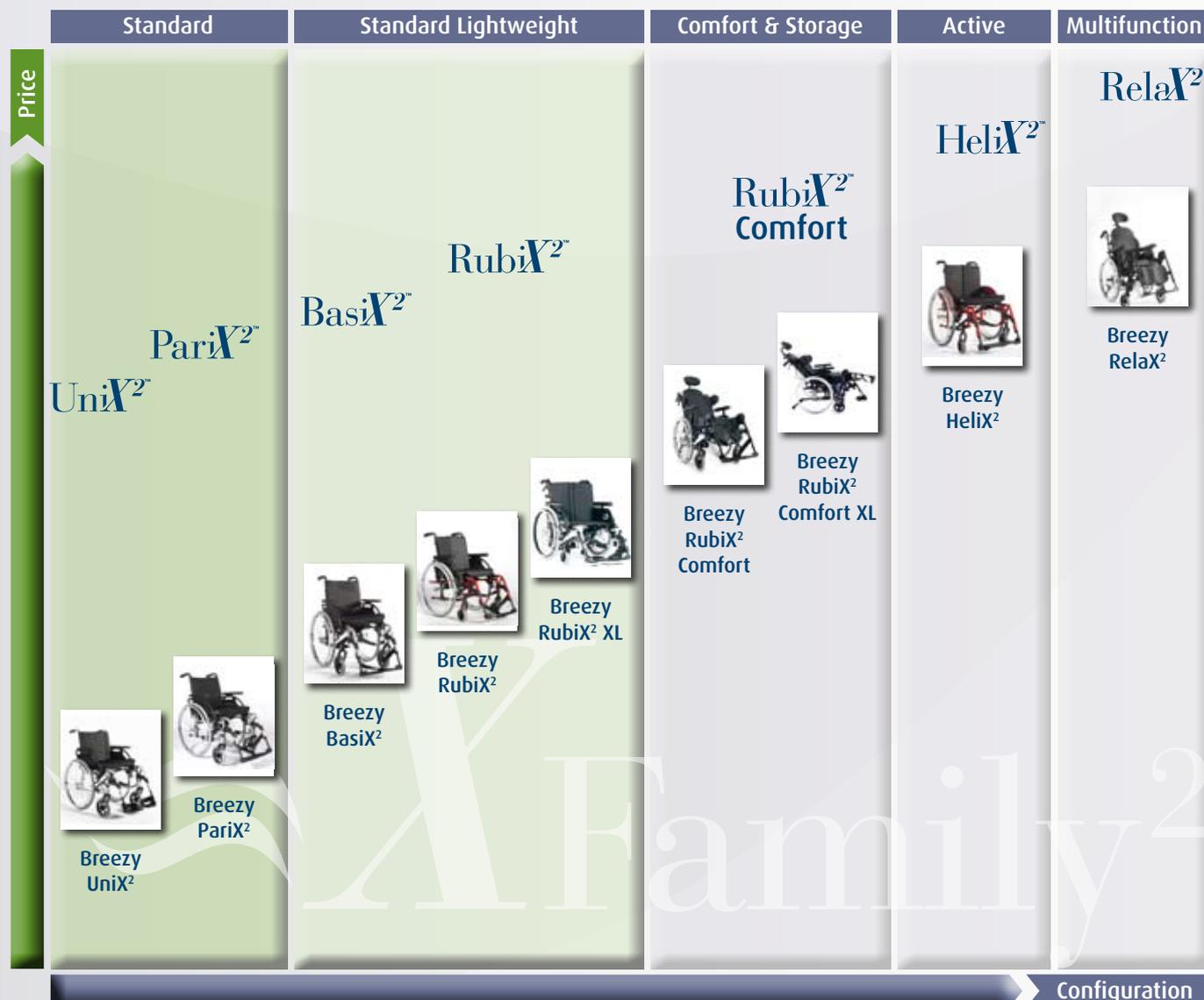

# Maximum use of common parts for optimised economy and perfect fleet management.

Surprisingly simple and yet so ingenious: The Breezy XFamily<sup>2</sup> Mini-Max principle. The high level of common parts means that you need to hold less stock while still offering the best in flexibility.

	Breezy UniX <sup>2</sup>	Breezy PariX <sup>2</sup>	Breezy BasiX <sup>2</sup>	Breezy RubiX <sup>2</sup>
Footrest	●	●	●	●
Elevating legrest	●	●	●	●
Stump pad	●	●	●	●
Armrest	-	●	●	●
Height adj. sideguard Standard version	●	●	●	●
Height adj. sideguard Central support	-	-	●	●
Wheel lock	●	●	●	●
Brake lever extension	●	●	●	●
6"/8" castors	●	●	●	●
Std. 22" / 24" wheels	●	●	●	●
22"/24" wheels drum brake	●	●	●	●
HeliX seat sling	●	-	-	-
BasiX seat sling	-	●	●	●
Anti-tip tubes	●	●	●	●
Spoke guard	●	●	●	●
Positioning belt	●	●	●	●

- Minimal handling effort, maximum interchangeability.
- Minimal stock levels for replacement parts, maximum individuality.
- Minimal order formalities, maximum profitability.

The Breezy XFamily<sup>2</sup> Mini-Max principle will convince you not only by the range available and the design, but also because it pays off for you too.

You stock your fleet according to your individual requirements, then we make sure that any subsequent costs over the whole product life cycle remain low.

	Breezy RubiX <sup>2</sup> XL	Breezy RubiX <sup>2</sup> Comfort	Breezy RubiX <sup>2</sup> Comfort XL	Breezy Helix <sup>2</sup>	Breezy RelaX <sup>2</sup>
	●	●	●	●	●
	●	●	●	●	●
	●	●	●	●	●
	●	●	●	●	-
	●	●	●	●	-
	●	●	●	●	●
	●	●	●	●	-
	●	●	●	●	-
	●	●	●	●	●
	●	●	●	●	●
	●	●	●	●	●
	-	-	-	-	-
	●	●	●	●	-
	●	●	●	●	●
	●	●	●	●	●
	●	●	●	●	●

# XFamily<sup>2</sup>





# UniX<sup>2</sup>

Build to last. Robust reliable design with full guarantee. Ensuring mobility both inside and outside.



Flip-up sideguard with adjustable arm-pad long and short.



Two adjustable seat heights.



Elevating legrest (option).

XFamily<sup>2</sup>

PariX<sup>2</sup><sup>™</sup>



Combination of steel and aluminium for higher rigidity.  
Widest adjustability in class.



Axle plate enables flexibility for seat heights and rear wheel position.



Adjustable castor angle ( 0° & 3°).  
Outrigger castor attachment.



Short and long comfortable armpad in one.



# BasiX<sup>2</sup><sup>™</sup>

Improved adjustability and rigidity combined with a more modern look!



New improved frame offers a better look and more seat depth adjustability (41 to 46 cm).



Tension adjustable backrest now extendable from 41 to 46 cm.



New robust castor attachment offers more rigidity and a better look ( 0°, 3° & 6°).

# RubiX<sup>2</sup><sup>™</sup>

Broad range of clinical options including hemiplegic specification, angle adjustable & reclining backs, comfort seating and Jay allows the chair to match needs of a wide user range.



New modern looking and robust 3 arm cross brace.



Comfort seat and reclining back (option) allows long and comfortable seating.



Tension adjustable 4 strap backrest extendable from 41 to 46 cm.



# RubiX<sup>2</sup> XL

The RubiX<sup>2</sup> XL extends the user range for heavy duty users by offering seat width up to 60 cm with 170 kg user weight.



New robust castor attachment offers more rigidity and a better look ( 0°, 3° & 6°).



Robust reliable 4 arm cross brace.



Drum brake (option).

# HeliX<sup>2</sup><sup>™</sup>



The configurable active chair designed for efficient fleet use which combines the strength of the RubiX<sup>2</sup> with lifestyle options.



Tension adjustable back sling and back height adjustment from 35-50 cm for optimal active ergonomics and comfort (standard feature).



Angle adjustable back to match active and semi-active user requirements (option).



Choice of lightweight stylish wheels, castors and tyres for the active user.



# Relax<sup>2</sup>

Modern looks, high comfort and functionality. A wheelchair which offers with its flexibility, adjustability and wide option range a perfect mobility solution .



Infinitely reclining backrest for optimum positioning. Recline from 0° - 35°. Seat tilt from -3° - 25°.



Infinitely adjustable seat depth to accommodate developing needs. Adjustable from 42 - 50 cm.



New and improved ELR for better functionality and look.

## Accessories and seating

The XFamily<sup>2</sup> comes standard with a wide range of options, settings and adjustments.



Height adjustable push handles



Reclining back (30° infinite adjustable)



Half folding back



Anti tips (swing away ; adjustable)




Drum brake



Height adjustable armrest with sliding pad for long and short position



Elevating legrest

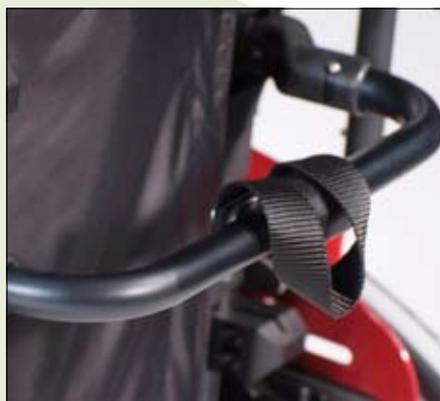

Foldable stabilizer bar



Stump pad



The XFamily<sup>2</sup> is designed so that the user can be provided with optimal comfort with different seats and backs. The comfort seats and backs are excellent to suit all individual needs. With the Jay range it is possible to optimize the user to position him properly and to support for his posture.



The Jay Basic is designed for the symmetrical client with low skin risk who requires comfort and mild postural support.



The Jay Soft Combi P cushions is designed for the client who has symmetrical posture and low risk of skin breakdown needing minimum to moderate postural support.



The Jay Easy Visco is a lightweight cushion with high level comfort and stability.



The Jay 3 Back is anthropometrically designed for support and positioning. Wide range of shapes and sizes to match user's clinical needs.

# The Breezy XFamily²!

Maximum use of common parts,  
optimised economy.



	Breezy UniX <sup>2</sup>	Breezy PariX <sup>2</sup>
Seat width	38 - 50 cm	38 - 52 cm
Seat depth	42 cm	42.5 cm
Seat height	47 - 51 cm	46 - 51 cm
Seat angle	3°	0° and 3°
Back height	40 - 42 cm	41 - 43 cm
Back angle	12°	12°
Safe slope	10°	10°
Centre of gravity adjustment	-	3,5 cm to -3,5 cm
Adjustable back upholstery	-	Yes
Sideguard	short and long	Tool height adjustable, short and long
Colour	Brilliant silver	Brilliant silver
Weight	from 18 kg	from 16.6 kg
Max. user weight	125 kg	125 kg



	<b>Breezy Basix<sup>2</sup></b>	<b>Breezy Rubix<sup>2</sup></b>	<b>Breezy Rubix<sup>2</sup> XL</b>
	38 - 52 cm	38 - 52 cm	52 - 60 cm
	43,5 cm (41 - 46 cm adjustable)	41 - 51 cm	41 - 51 cm
	42,5 - 51 cm	37,5 - 53 cm	37,5 - 51 cm
	0° to 6°	0° to 6°	0° to 6°
	43,5 cm (41 - 46 cm adjustable)	43,5 cm (41 - 46 cm adjustable)	43,5 cm (41 - 46 cm adjustable)
	12° (7° - 30° with reclining back)	12° (-5° - 20° with angle adj back)	12° (-5° - 20° with angle adj back)
	10°	10°	10°
	5 cm to -5 cm	5 cm to -5 cm	5 cm to -5 cm
	Yes	Yes	Yes
	Tool height adjustable, short and long	Tool height adj., short and long sliding	Tool height adj., short and long sliding
	Brilliant silver	Brilliant silver, blue, brilliant red	Brilliant silver, blue, brilliant red
	from 14,9 kg	from 15,5 kg	from 16,8 kg
	125 kg	125 kg	170 kg

# The Breezy XFamily²!

Maximum use of common parts,  
optimised economy.



	<b>Breezy RubiX² Comfort</b>	<b>Breezy Helix²</b>
Seat width	38 - 52 cm (52 - 60 cm XL-Version)	38 - 52 cm (52 - 60 cm XL-Version)
Seat depth	41 - 48,5 cm	41 - 51 cm
Seat height	37,5 - 53 cm (37,5 - 51 cm XL-Version)	35.5 - 53 cm (37,5 - 51 cm XL-Version)
Seat angle	0° to 6°	0° to 6°
Back height	46 cm	35 - 50 cm
Back angle	7° to 30° adjustable	12° (-5° to 20° angle adjustable)
Safe slope	10°	10°
Centre of gravity adjustment	5 cm to -5 cm	-5 cm to 7.5 cm
Adjustable back upholstery	Yes (Comfort back optional)	Yes
Sideguard	Tool height adj., short and long, sliding	height adjustable, short and long
Colour	Brilliant silver, blue, brilliant red	Brilliant silver, blue, brilliant red, black
Weight	from 15,5 kg / XL Version from 16,8 kg	from 15.3 kg / XL Version from 16.3 kg
Max. user weight	125 kg / XL-Version 170 kg	125 kg / XL-Version from 170 kg



**Breezy  
RelaX<sup>2</sup>**

41(37 - 43), 46(42 - 48), 51(47 - 53) cm, 53 cm XL-Version)

42 - 50 cm (infinite)

42.5 - 50 cm

0° to 25° / -3°- 22°

55 - 60 cm

0° to 35° (infinite)

10°

-


-

height- /depth adjustable, sliding

Brilliant silver

from 34 kg

145 kg (170 kg XL-Version)

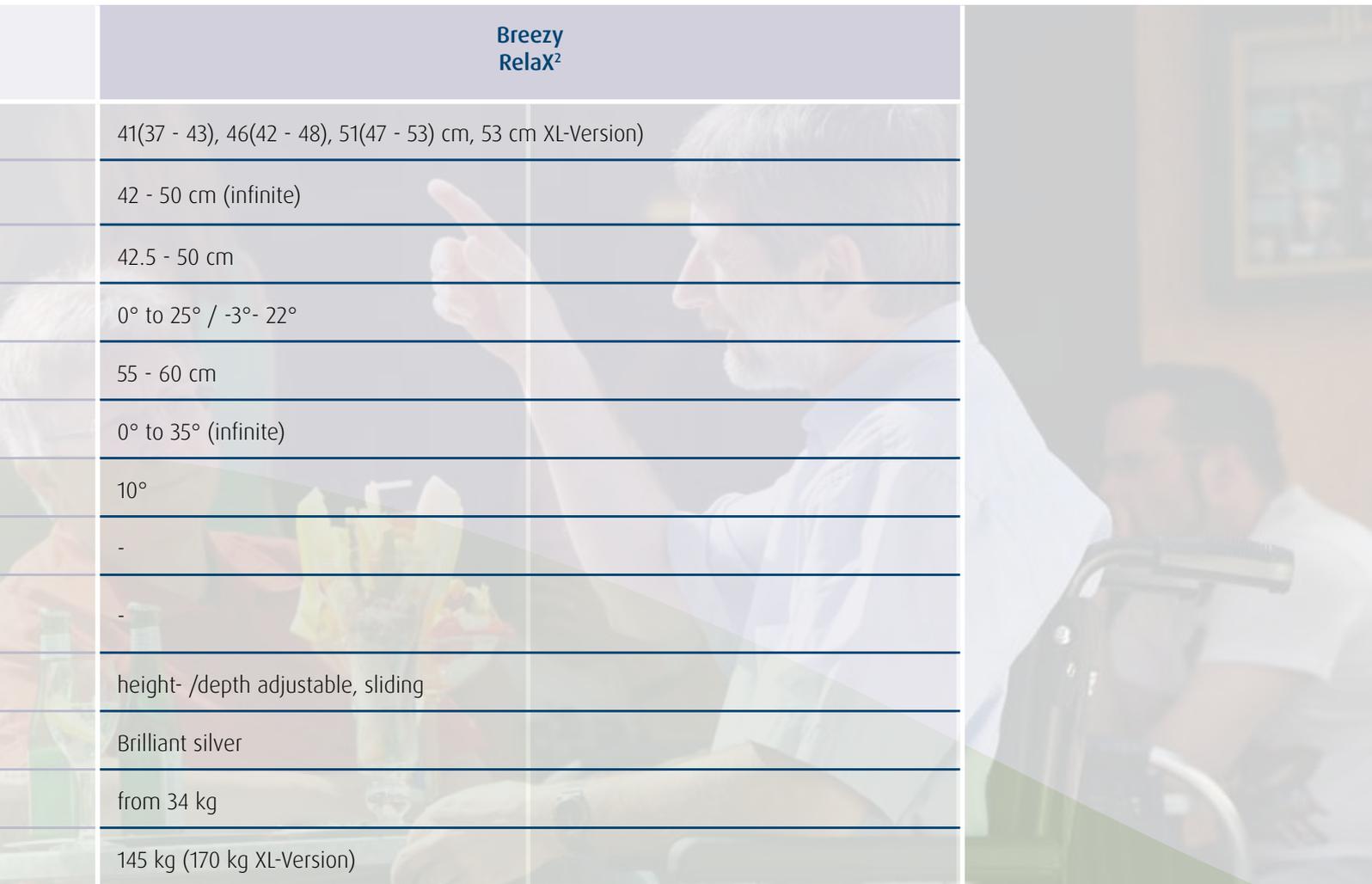



 **Breezy**<sup>®</sup>  
Premium products. Quality service.

All products in the Breezy XFamily<sup>2</sup> are crash-tested in accordance with ISO 7176-19 and are therefore approved for the transportation of a person in a vehicle, whilst seated in the wheelchair.



**Sunrise Medical GmbH**

Kahlbachring 2-4  
69254 Malsch/HD  
Tel.: +49 (0) 72 53 / 9 80-400  
Fax: +49 (0) 72 53 / 9 80-222

E-mail: [international@sunrisemedical.de](mailto:international@sunrisemedical.de)  
[www.SunriseMedical.eu](http://www.SunriseMedical.eu)